



FreeStand_™



Enhancing Outcomes for Patients and Their Caregivers...

A freestanding alternative.





FreeStand_{TM} is a new versatile freestanding lift system, well worth its price. Easy to install and fits most environments. Since FreeStand_{TM} is completely freestanding, there are no special requirements on the ceiling or on the walls.

The reasons to choose a freestanding system may vary: for example the system is to be tried only for a limited period of time or the lifting need is temporary. Sometimes a fixed installation is not possible to perform because of the structure of the room.

FreeStand's leg supports and the rail are white, to blend into most private homes and care facilities. The legs have a wide, stable base and the system has smart support arms that can be streched out against the wall to provide extra stability.

The rail is available in several different lengths between 2.0 and 3.2 metres (6' 6" to 10' 6").

For a complete system:

3104001 FreeStand_™ leg supports 3104002 Parking unit

Rail H100 up to 3,2 metres (10' 6")

Lift motor: Likorall™ or Multirall™

Carriage Hand Control Charger

Lifting Accessory




Support arms can be stretched against the wall to provide extra stability.



A complete system consists of two leg supports (sold as one article), a rail, a lift motor with carriage, hand control and charger, and a suitable lifting accessory.



The charger can be mounted on one of the leg supports. The parking position for the slingbar is also easily mounted here.

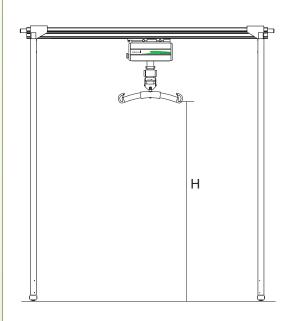

Lifting height with lift motor


Maximum lifting height (H) depends on which lift motor that is used.


Examples:

- Likorall_{**} 242 with mechanical emergency lowering: H = 1715 mm (5' 7")
- Likorall[∞] 242 without mechanical emergency lowering: H = 1820 mm (5' 11")
- Multirall∞: H = 1740 mm (5' 8") All examples with Universal Slingbar 450.

With a common sling (UniversalSling™, mod 00 or Original HighBack Sling™, mod 20), the total lifting height (from the floor to the lower edge of the sling) is appr. 70-80 centimeters above the floor (2′ 5″).

Contact Liko for more information.

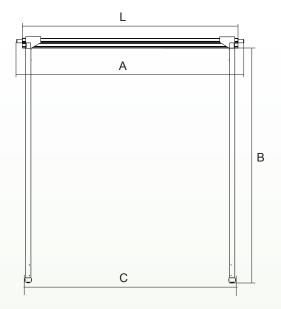

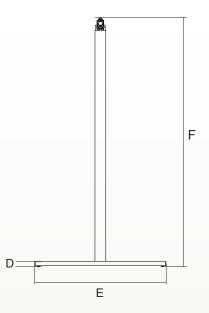

Max. load¹	L	A max	A min	В	С	D	Е	F	Weight	
									tot	part*
200/400	2000/79	3710/146	2040/81	2180/85	2010/80	45/2	121/5	2302/91	27,5/60.6	7,3/16.1
200/400	2200/87	3910/153	2240/89	2180/85	2210/87	45/2	121/5	2302/91	28,2/62.2	8,0/17.6
200/400	2400/94	4110/161	2440/97	2180/85	2410/95	45/2	121/5	2302/91	29,0/63.9	8.8/19.4
200/400	2600/102	4310/169	2640/104	2180/85	2610/103	45/2	121/5	2302/91	29,7/65.5	9,5/20.9
200/400	2800/110	4510/177	2840/112	2180/85	2810/111	45/2	121/5	2302/91	30,4/67.0	10,2/22.5
200/400	3000/118	4710/185	3040/120	2180/85	3010/119	45/2	121/5	2302/91	31,2/68.8	11,0/24.3
200/400	3200/126	4910/193	3240/128	2180/85	3210/127	45/2	121/5	2302/91	31,9/70.3	11,7/25.8

Max. load and weight in kg/lbs. Dimensions in mm/inch.

* Weight of the single heaviest component on the FreeStand™ system when disassembled.

((



FreeStand $_{100}$ complies with applicable requirements of the Medical Device Directive for Class-1 products (MDD 93/42/EEC). FreeStand $_{100}$ complies with EN ISO 10535.

Max. load: 200 kg (440 lbs.)

Design and Quality by Liko in Sweden.

Liko is quality-certified in accordance with ISO 9001 and its equivalent for the medical device industry, ISO 13485.

Liko is also certified in accordance with environmental standard ISO 14001.



Enhancing Outcomes for Patients and Their Caregivers.™

France +33 (0)2 97 50 92 12
United Kingdom +44 (0)1530 411000
Deutschland +49 (0)211 16450-0
Nederland +31 (0)347 32 35 32
Italia +39 (0)2 950541
Suisse/Schweiz +41 (0)21 706 21 30
(deutschsprachig) +41 (0)21 706 21 38

 Österreich
 +43 (0)2243 28550

 Ireland
 +353 (0)1 413 6005

 Iberia
 +34 (0)93 6856000

 Nordic Countries
 +46 (0)20 78 10 30

 Export
 Contact your local Hill-Rom distributor or contact your Area Manager via website or call
 +1 812 934 8173

www.hill-rom.com www.liko.com

Hill-Rom_® reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom_® makes is the express written warranty extended on the sale or rental of its products. © 2012 Hill-Rom_® SARL, ALL RIGHTS RESERVED.

Doc. No. 7EN131301-02 • Design: Hill-Rom International Communications, 2012-11